

SunFire Columns

CONTENTS

I. GETTING STARTED

- a. Column Installation
- b. Column Equilibration
- c. Initial Column Efficiency Determination

II. COLUMN USE

- a. Guard Columns
- b. Sample Preparation
- c. pH Range
- d. Solvents
- e. Pressure
- f. Temperature

III. SCALING UP/DOWN ISOCRATIC METHODS

IV. TROUBLESHOOTING

V. COLUMN CLEANING, REGENERATING AND STORAGE

- a. Cleaning and Regenerating
- b. Storage

VI. CONNECTING THE COLUMN TO THE HPLC SYSTEM

- a. Column Connectors and System Tubing Considerations
- b. Bandspreading Minimization
- c. Measuring System Bandspreading Volume and System Variance
- d. Measuring System Volume

VII. ADDITIONAL INFORMATION

- a. Availability and Use of Electronic Column Tags
- b. Use of Narrow-Bore (≤ 3.0 mm i.d.) Columns
- c. Impact of Bandspreading Volume on 2.1 mm i.d. Column Performance
- d. Non-Optimized vs. Optimized LC/MS/MS System: System Modification Recommendations

Thank you for choosing a SunFire™ Column. The SunFire packing materials were designed to provide excellent peak shape, minimal column bleed, and high mass loading. The SunFire packing materials are manufactured in an ISO 9000 certified plant using ultra pure reagents. Each batch of SunFire material is tested chromatographically with acidic, basic and neutral analytes and the results are held to narrow specification ranges to assure excellent, reproducible performance. Every column is individually tested and a Performance Chromatogram is provided with each column along with the Certificate of Batch Analysis.

SunFire Column Physical Characteristics

Chemistry	C ₁₈	C ₈	Silica
USP class no.	L1	L7	L3
Particle shape	Spherical	Spherical	Spherical
Particle sizes	2.5, 3.5, 5, 10 μm	2.5, 3.5, 5, 10 μm	5, 10 μm
Pore size	100 Å	100 Å	100 Å
Surface area	340 m ² /g	340 m ² /g	340 m ² /g
Carbon load	16%	12%	N/A
Endcapped	Yes	Yes	N/A
pH range	2–8	2–8	N/A
Temperature limits	Low pH = 50 °C High pH = 40 °C	Low pH = 40 °C High pH = 40 °C	N/A



I. GETTING STARTED

Each SunFire Column comes with a Performance Test Chromatogram. This Performance Test Chromatogram is specific to each individual column and contains the following information: gel batch number, column serial number, USP plate count, USP tailing factor, capacity factor, and chromatographic conditions. The performance test chromatogram should be stored for future reference.

a. Column Installation

Note: Flow rates given in the procedure below are for a typical 4.6 mm i.d. column. Scale the flow rate up or down accordingly based upon the column i.d., length, particle size, and backpressure of the SunFire Column being installed. See "Scaling Up/Down Isocratic Separations" for calculating flow rates when changing column i.d. and/or length. See "Connecting the Column to the HPLC System" for a more detailed discussion on HPLC connections.

Reversed-phase Columns (SunFire C₁₈ and SunFire C₈)

1. Purge the pumping system of any buffer-containing mobile phases and connect the inlet end of the column to the injector outlet.
2. Flush column with 100% organic mobile phase (methanol or acetonitrile) by setting the pump flow rate to 0.1 mL/min and increase the flow rate to 1 mL/min over 5 minutes.
3. When the mobile phase is flowing freely from the column outlet, stop the flow and attach the column outlet to the detector. This prevents entry of air into the detection system and gives more rapid baseline equilibration.
4. Gradually increase the flow rate as described in step 2.
5. Once a steady backpressure and baseline have been achieved, proceed to the next section.

Normal-phase Columns (SunFire Silica)

Note: It is assumed that your system has been used for reversed-phase chromatography. If this is not the case, you can start with step 3.

1. Purge the pumping system of any buffer-containing mobile phases.
2. Flush the system thoroughly with acetonitrile.
3. Switch the system over to the mobile phase that you are planning to use in normal-phase chromatography.
4. Connect the column and equilibrate it with the mobile phase.

Note: Equilibration with the mobile phase may require a larger amount of solvent than in reversed-phase chromatography.

b. Column Equilibration

SunFire columns are shipped in 100% acetonitrile.

It is important to ensure mobile phase compatibility before changing to a different mobile-phase system. Equilibrate the column with a minimum of 10 column volumes of the mobile phase to be used (refer to Table 1 for a listing of empty column volumes).

Reversed-phase (SunFire C₁₈ or SunFire C₈) Columns

To avoid precipitating out mobile-phase buffers on your column or in your system, flush the column with five column volumes of a water/organic solvent mixture, using the same or lower solvent content as in the desired buffered mobile phase. (For example, flush the column and HPLC system with 60% methanol in water prior to introducing 60% methanol/40% buffer mobile phase.)

Note: If mobile-phase additives are present in low concentrations (e.g., ion-pairing reagents), 100 to 200 column volumes may be required for complete equilibration.

Normal-phase Columns (SunFire Silica)

SunFire normal-phase (NP) columns are delivered in 96% heptane/4% isopropyl alcohol. Care should be taken not to pass any mobile phase through the column that might cause a precipitate (see above). SunFire Silica normal-phase columns are compatible with water and all common organic solvents, provided that solvent miscibility is accounted for.

Equilibrate normal-phase silica columns in the mobile phase. Very small quantities of water in the mobile phase can dramatically affect the activity of normal-phase packings. For good reproducibility, ensure that the mobile phase always has the same water content.

It is difficult and usually unnecessary to completely eliminate the water from the mobile phase. Dry mobile phases can take a very long time to equilibrate the column. A water content of 50 percent of saturation is recommended for most applications.

To equilibrate your column:

1. Starting at 0.0 mL/min, increase the flow rate in 0.1 mL/min increments to 1.0 min.
2. Purge the column with the mobile phase until you obtain a stable baseline.
3. Verify that retention times and peak areas for a standard are stable by comparing 2–3 replicate, consecutive injections.

Before you perform the first analysis on your new column, perform an efficiency test to confirm the performance of the column.

c. Initial Column Efficiency Determination

1. Perform an efficiency test on the column before using it. Waters recommends using a suitable solute mixture, as found in the "Performance Test Chromatogram", to analyze the column upon receipt. However, if the column is used only for a single routine isocratic assay, it may be more convenient to test the column under these assay conditions. Keep a record of the initial column performance.
2. Determine the number of theoretical plates (N) and use this value for periodic comparisons.
3. Repeat the test at predetermined intervals to track column performance over time. Slight variations may be obtained on two different HPLC systems due to the quality of the connections, operating environment, system electronics, reagent quality, column condition and operator technique.

Table 1. Empty Column Volumes in mL (multiply by 10 for flush solvent volumes)

Column Length	Column Internal Diameter (nm)								
	1.0	2.1	3.0	4.6	7.8	10	19	30	50
20 mm	-	0.07	0.14	0.33	-	-	-	-	-
30 mm	-	0.1	0.2	0.5	-	2.4	8	-	-
50 mm	0.1	0.2	0.3	0.8	2.4	4	14	35	98
100 mm	0.1	0.4	0.7	1.7	5	8	28	70	196
150 mm	0.1	0.5	1.0	2.5	7	12	42	106	294
250 mm	-	0.9	1.8	4	-	20	70	176	490

II. COLUMN USE

To ensure the continued high performance of SunFire Columns, follow these guidelines:

a. Guard Columns

Use a SunFire guard column of matching chemistry and particle size between the injector and main column. It is important to use a matching guard column to protect the main column while not compromising or changing the analytical resolution. Guard columns need to be replaced at regular intervals as determined by sample contamination. When system backpressure steadily increases above a set pressure limit, it is usually an indication that the guard column should be replaced. A sudden appearance of split peaks or other changes in chromatographic performance is also indicative of a need to replace the guard column.

b. Sample Preparation

1. Sample impurities often contribute to column contamination. One option to avoid this is to use Oasis™ Solid-phase Extraction Cartridges/Columns or Sep-Pak™ Cartridges of the appropriate chemistry to clean up the sample before analysis. Link to waters.com/sampleprep
2. It is preferable to prepare the sample in the operating mobile phase or a mobile phase that is weaker (less organic modifier in the case of reversed-phase chromatography, less polar modifier in the case of normal-phase chromatography or hydrophilic interaction chromatography, less salt in the case of ion exchange) than the mobile phase for the best peak shape and sensitivity.
3. If the sample is not dissolved in the mobile phase, ensure that the sample, solvent and mobile phases are miscible in order to avoid sample and/or buffer precipitation. Filter sample with 0.2 µm membranes to remove particulates. If the sample is dissolved in a solvent that contains an organic modifier (e.g., acetonitrile, methanol, etc.) ensure that the membrane material does not dissolve in the solvent. Contact the membrane manufacturer with solvent compatibility questions. Alternatively, centrifugation for 20 minutes at 8,000 rpm, followed by the transfer of the supernatant liquid to an appropriate vial, could be considered.

c. pH Range

The recommended operating pH range for SunFire Columns is 2 to 8. A listing of commonly used buffers and additives is given in Table 2. Additionally, the column lifetime will vary depending upon the operating temperature, and the type and concentration of buffer used. For example, the use of phosphate buffer at pH 8 in combination with elevated temperatures will lead to shorter column lifetimes.

Table 2: Buffer Recommendations for Using SunFire Columns from pH 2 to 8

Additive or Buffer	pK _a	Buffer Range (±1 pH unit)	Volatility	Used for Mass Spec?	Comments
Formic Acid	3.75		Volatile	Yes	Maximum buffering obtained when used with Ammonium Formate salt. Used in 0.1–1.0% range.
Acetic Acid	4.76		Volatile	Yes	Maximum buffering obtained when used with Ammonium Acetate salt. Used in 0.1–1.0% range.
Ammonium Formate (NH ₄ COOH)	3.75	2.75–4.75	Volatile	Yes	Used in the 1–10 mM range for LC/MS. Higher concentrations (typically 20 mM) are recommended for UV applications. <i>Note: sodium or potassium salts are not volatile.</i>
Trifluoroacetic Acid (TFA)	0.30		Volatile	Low conc.	When used in LC/MS, due to signal suppression, it is generally recommended to use TFA at concentrations < 0.1%.
Ammonium Acetate (NH ₄ CH ₃ COOH)	4.76	2.75–5.76	Volatile	Yes	Used in the 1–10 mM range for LC/MS. Higher concentrations (typically 20 mM) are recommended for UV applications. <i>Note: sodium or potassium salts are not volatile.</i>
Phosphate 1	2.15	1.15–3.15	Non-volatile	No	Traditional low pH buffer, good UV transparency.
Phosphate 2	7.20	6.20–8.20	Non-volatile	No	Above pH 7, reduce temperature/concentration and use guard column to maximize lifetime.

d. Solvents

To maintain maximum column performance, use high quality chromatography grade solvents. Filter all aqueous buffers prior to use. Pall Gelman Laboratory Acrodisc® Filters are recommended. (Please refer to the filtration section of the Waters Chromatography Columns and Supplies Catalog or the Waters web site (waters.com) for additional product information.) Solvents containing suspended particulate materials will generally clog the outside surface of the inlet distribution frit of the column. This will result in higher operating pressure and poorer performance. Degas all solvents thoroughly before use to prevent bubble formation in the pump and detector. The use of an on-line degassing unit is also recommended. This is especially important when running low pressure gradients since bubble formation can occur as a result of aqueous and organic solvent mixing during the gradient.

e. Pressure

SunFire columns can tolerate pressures of up to 6,000 psi (400 bar or 40 MPa) although long-term, routine operating pressures greater than 4,000–5,000 psi should be avoided in order to maximize column and system lifetimes.

f. Temperature

Temperatures between 20 °C–50 °C are recommended for operating SunFire columns at low pH in order to enhance selectivity, lower solvent viscosity and increase mass transfer rates. However, any temperature above ambient may have a negative effect on lifetime which will vary depending on the pH and buffer conditions used.

III. SCALING UP/DOWN ISOCRATIC METHODS

The following formulas will allow scale up or scale down, while maintaining the same linear velocity, and provide new sample loading values:

If column *i.d.* and length are altered:

$$F_2 = F_1(r_2/r_1)^2$$

or

$$\text{Injection volume}_1 = \text{Injection volume}_2 (r_2/r_1)^2 (L_2/L_1)$$

Where: *r* = radius of the column, in mm

F = flow rate, in mL/min

L = length of column, in mm

1 = original, or reference column

2 = new column

IV. TROUBLESHOOTING

Changes in retention time, resolution, or backpressure are often due to column contamination. See "Column Cleaning, Regenerating and Storage". A copy of the HPLC Troubleshooting Guide may be downloaded at waters.com, in the "Search" field enter WA20769.

V. COLUMN CLEANING, REGENERATING, AND STORAGE

a. Cleaning and Regenerating

Changes in peak shape, peak splitting, shoulders on the peak, shifts in retention, change in resolution or increasing backpressure may indicate contamination of the column. Changing the guard column being used will often restore column performance. If not (or if no guard column is being used), follow the procedures detailed below. To prevent potential contamination from affecting detector performance, it is recommended that any detector(s) be disconnected from the effluent flow of the column during cleaning. Reversing the direction of the flow through the column (backflushing) may sometimes improve the effectiveness of any cleaning procedure.

Reversed-phase Columns (SunFire C₁₈ and SunFire C₈)

Flushing with a neat organic solvent, taking care not to precipitate buffers, is usually sufficient to remove most contaminant. If the flushing procedure does not solve the problem, purge the column using the following cleaning and regeneration procedures. Use the cleaning routine that matches the properties of the samples and/or what you believe is contaminating the column (see Table 3). Flush columns with 20 column volumes of HPLC-grade solvents (e.g., 80 mL total for 4.6 x 250 mm column). Increasing mobile-phase temperature to 35–55 °C increases cleaning efficiency. If the column performance is poor after regenerating and cleaning, call your local Waters office for additional support.

Table 3: Column Sequence or Options

Polar Samples	Non-polar Samples	Proteinaceous Samples
1. Water	1. Isopropanol (or an appropriate isopropanol/water mixture*)	Option 1: Inject repeated aliquots of dimethyl sulfoxide (DMSO)
2. Methanol	2. Tetrahydrofuran (THF)	Option 2: Gradient of 10% to 90% B where:
3. Tetrahydrofuran (THF)	3. Dichloromethane	A = 0.1% trifluoroacetic acid (TFA) in Water
4. Methanol	4. Hexane	B = 0.1% trifluoroacetic acid (TFA) in acetonitrile (CH ₃ CN)
5. Water	5. Isopropanol (followed by an appropriate isopropanol/water mixture*)	Option 3: Flush column with 7 M guanidine hydrochloride, or 7 M urea
6. Mobile Phase	6. Mobile Phase	

*Use low organic solvent content to avoid precipitating buffers.

Normal-phase Columns (SunFire Silica)

To regenerate, pump 20–30 column volumes each of dichloromethane and isopropanol through the column. Other wash solvents such as tetrahydrofurane (THF) may also be selected based on the suspected contamination.

Guard columns need to be replaced at regular intervals, as determined by sample contamination. When system backpressure steadily increases above a set pressure limit, it is usually an indication that the guard column should be replaced. A sudden appearance of split peaks is also indicative of a need to replace the guard column.

b. Storage

Completely seal the column to avoid evaporation and drying out of the bed.

Reversed-phase Columns (SunFire C₁₈ and SunFire C₈)

For periods longer than four days at room temperature, store the column (with the exception of cyano chemistry columns) in 100% acetonitrile at room temperature. For elevated temperature applications, store immediately after use in 100% acetonitrile for the best column lifetime. Do not store columns in buffered eluents. If the mobile phase contained a buffer salt, flush the column with 10 column volumes of HPLC grade water (see Table 1 for common column volumes) and replace with 100% acetonitrile for storage. Failure to perform this intermediate step could result in precipitation of the buffer salt in the column when 100% acetonitrile is introduced.

Normal-phase Columns (SunFire Silica)

For rapid equilibration upon start-up, store your normal-phase column in the mobile phase that is commonly used.

VI. CONNECTING THE COLUMN TO THE HPLC SYSTEM

a. Column Connectors and System Tubing Considerations

All SunFire Columns have Waters-style endfittings.

Tools needed for SunFire analytical column:

- 1/2 inch wrench
- 9/16 inch wrench

Handle the column with care. Do not drop or hit the column on a hard surface as it may disturb the bed and affect its performance.

1. Correct connection of 1/16 inch outer diameter stainless steel tubing leading to and from the column is essential for high quality chromatographic results.
2. When using standard stainless steel compression screw fittings, it is important to ensure proper fit of the 1/16 inch outer diameter stainless steel tubing. When tightening or loosening the compression screw, place a 5/16 inch wrench on the compression screw and a 1/2 inch wrench on the hex head of the column endfitting.

Note: If one of the wrenches is improperly placed on the column tube flat during this process, the endfitting will be loosened and leak.

3. If a leak occurs between the stainless steel compression screw fitting and the column endfitting, a new compression screw fitting, tubing and ferrule must be assembled.
4. An arrow on the column identification label indicates correct direction of solvent flow. Correct connection of 1/16 inch outer diameter stainless steel tubing leading to and from the column is essential for high-quality chromatographic results. It is important to realize that extra column peak broadening can destroy a successful separation. The choice of appropriate column connectors and system tubing is discussed in detail below. Due to the absence of an industry standard, various column manufacturers have employed different types of chromatographic column connectors. The chromatographic performance of the separation can be negatively affected if the style of the column endfittings does not match the existing tubing ferrule settings. This section explains the differences between Waters style and Parker style ferrules and endfittings (Figure 1). Each endfitting style varies in the required length of the tubing protruding from the ferrule. The SunFire column is equipped with Waters style endfittings that require a 0.130 inch ferrule. If a Parker style column is presently being used, it is critical that the ferrule depth be reset for optimal performance prior to installing a SunFire Column.

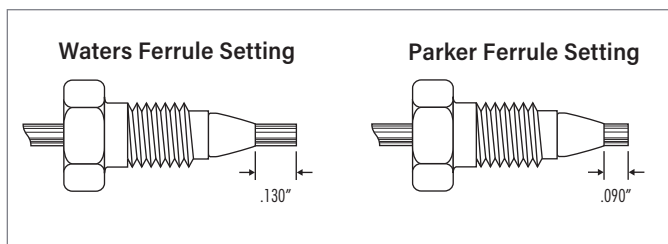


Figure 1: Waters and Parker Style Ferrule Types.

In a proper tubing/column connection (Figure 2), the tubing touches the bottom of the column endfitting, with no void between them.

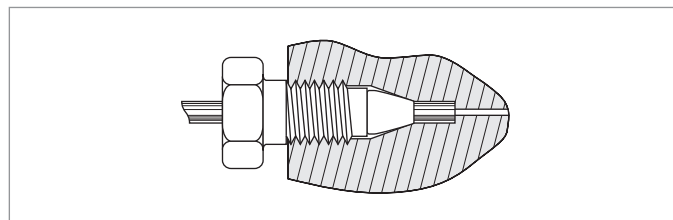


Figure 2: Proper Tubing/Column Connection.

The presence of a void in the flow stream reduces column performance. This can occur if a Parker style ferrule is connected to a Waters endfitting (Figure 3).

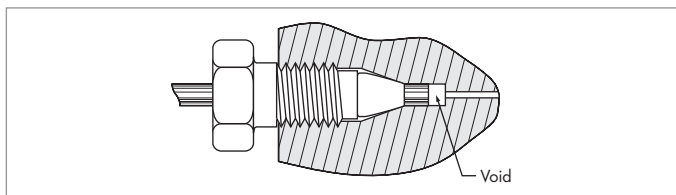


Figure 3: Parker Ferrule in a Waters Style Endfitting.

Note: A void appears if tubing with a Parker style ferrule is connected to a Waters style column.

There is only one way to fix this problem: Cut the end of the tubing with the ferrule, place a new ferrule on the tubing and make a new connection. Before tightening the screw, make sure that the tubing bottoms out in the endfitting of the column. Conversely, if tubing with a Waters ferrule is connected to a column with Parker style endfitting, the end of the tubing will bottom out before the ferrule reaches its proper sealing position. This will leave a gap and create a leak (Figure 4).

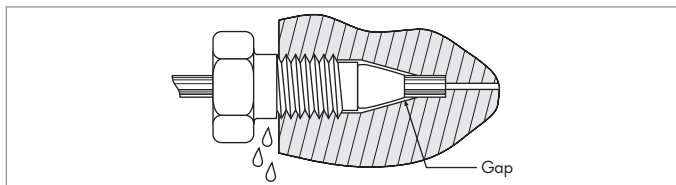


Figure 4: Waters Ferrule in a Parker Style Endfitting.

Note: The connection leaks if a Waters ferrule is connected to a column with a Parker style endfitting.

There are two ways to fix the problem:

1. Tighten the screw a bit more. The ferrule moves forward, and reaches the sealing surface. Do not overtighten since this may end in breaking the screw.
2. Cut the tubing, replace the ferrule and make a new connection. Alternatively, replace the conventional compression screw fitting with an all-in-one PEEK™ Fitting (Waters Part Number PSL613315) that allows resetting of the ferrule depth. (Note: PEEK Fittings are not recommended for normal-phase applications!) Another approach is to use a Thermo Corporation SLIPFREE™ Connector to always ensure the correct fit.

The fingertight SLIPFREE Connectors automatically adjust to fit all compression screw type fittings without the use of tools (Figure 5).

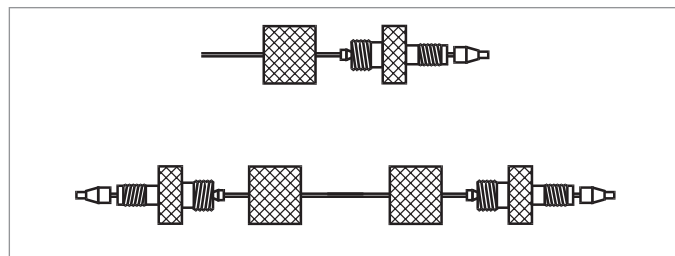


Figure 5: Single and Double SLIPFREE Connectors

SLIPFREE Connectors Features

- Tubing pushed into endfitting, thereby guaranteeing a void-free connection
- Connector(s) come(s) installed on tubing
- Various tubing i.d. and lengths available
- Fingertight to 10,000 psi – never needs wrenches
- Readjusts to all column endfittings
- Compatible with all commercially available endfittings
- Unique design separates tube-holding function from sealing function

Table 4. Waters Part Numbers for SLIPFREE Connectors

SLIPFREE® Type	Tubing Internal Diameter		
Tubing Length	0.005"	0.010"	0.020"
Single 6 cm	PSL 618000	PSL 618006	PSL 618012
Single 10 cm	PSL 618002	PSL 618008	PSL 618014
Single 20 cm	PSL 618004	PSL 618010	PSL 618016
Double 6 cm	PSL 618001	PSL 618007	PSL 618013
Double 10 cm	PSL 618003	PSL 618009	PSL 618015
Double 20 cm	PSL 618005	PSL 618001	PSL 618017

b. Bandspreading Minimization

Figure 6 shows the influence of tubing internal diameter on system bandspreading and peak shape. As can be seen, the larger tubing diameter causes excessive peak broadening and lower sensitivity.

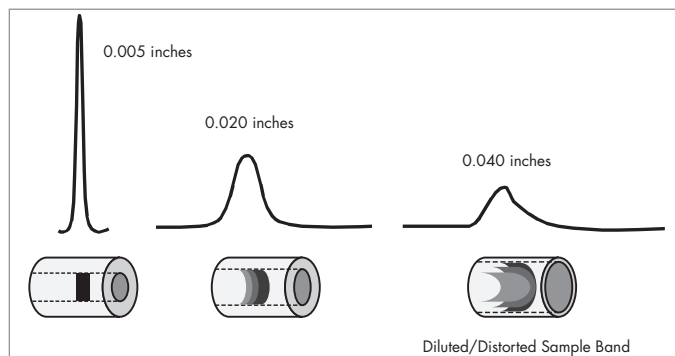


Figure 6: Effect of Connecting Tubing on System.

c. Measuring System Bandspreading Volume and System Variance

This test should be performed on an HPLC system with a single wavelength UV detector (not a Photodiode Array [PDA]).

1. Disconnect column from system and replace with a zero dead volume union.
2. Set flow rate to 1 mL/min.
3. Dilute a test mix in mobile phase to give a detector sensitivity of 0.5–1.0 AUFS (system start up test mix can be used which contains uracil, ethyl and propyl parabens; Waters Part Number WAT034544).
4. Inject 2 to 5 μL of this solution.
5. Measure the peak width at 4.4% of peak height (5-sigma method):

5-sigma Bandspreading (μL) = Peak Width (min) x Flow Rate (mL/min) x (1000 $\mu\text{L}/1\text{ mL}$)
 System Variance (μL^2) = (5-sigma bandspreading)² / 25

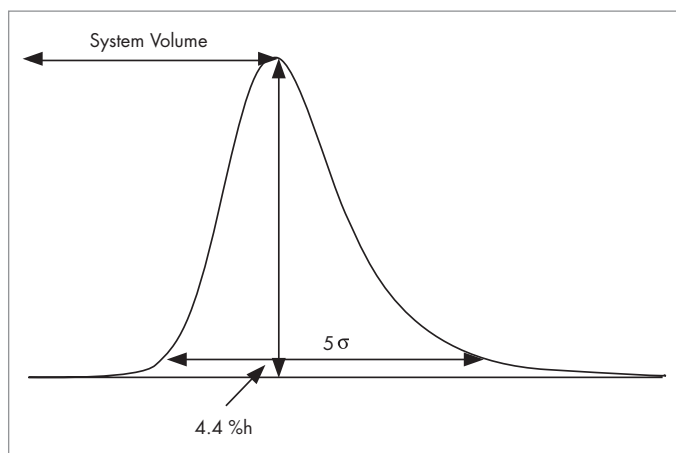


Figure 8: Determination of System Bandspreading Volume Using 5-Sigma Method.

In a typical HPLC system, the Bandspreading Volume should be $100\ \mu\text{L} \pm 30\ \mu\text{L}$ (or Variance of $400\ \mu\text{L}^2 \pm 36\ \mu\text{L}^2$).

In a microbore (2.1 mm i.d.) system, the Bandspreading Volume should be no greater than 20 to 40 μL (or Variance no greater than $16\ \mu\text{L}^2$ to $64\ \mu\text{L}^2$).

d. Measuring System Volume

System volume is important in scaling separations because it creates an isocratic hold at the start of every run. This hold is often several column volumes on a small scale, but a fraction of the volume of a prep column. Compensation for this volume must be included in planning a scaling experiment to avoid distorting the chromatography (Figure 9).

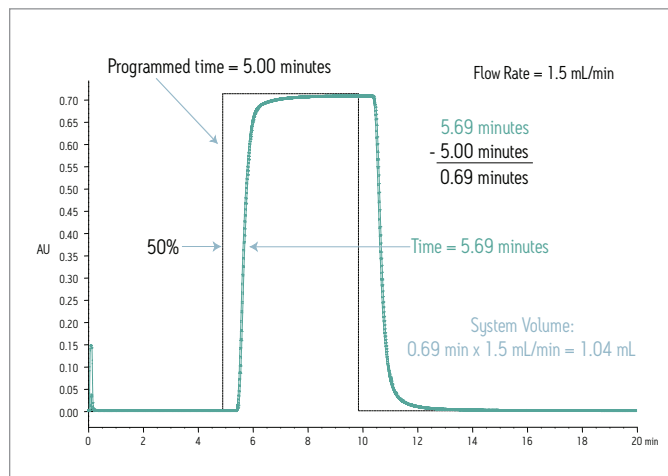


Figure 9: Determination of Gradient Delay Volume

1. Remove column.
2. Use acetonitrile as mobile phase A, and acetonitrile with 0.05 mg/mL uracil as mobile phase B (eliminates non-additive mixing and viscosity problems).
3. Set UV detector at 254 nm.
4. Use the flow rate in the original method and the intended flow rate on the target instrument.
5. Collect 100% A baseline for 5 minutes.
6. Program a step change at 5 minutes to 100% B, and collect data for an additional 5 minutes.
7. Measure absorbance difference between 100% A and 100% B.
8. Measure time at 50% of that absorbance difference.
9. Calculate time difference between start of step and 50% point.
10. Multiply time difference by flow rate.

VII. ADDITIONAL INFORMATION

a. Availability and Use of Electronic Column Tags

'Smart columns' are those enabled with electronic devices (tags) permanently affixed to enable identification and/or tracking of column history and use information. These tags contain unique product identifying information, such as part and serial numbers, and product descriptions. Data residing on these devices is automatically imported by complementary Waters' systems, whether LC instrumentation or the Empower™ Chromatographic Data System. This type of product data is used to aid in tracking and decision-making and troubleshooting, especially valuable for GLP/GMP compliant laboratories.

Tags come in two types: tethered like eCord™ and untethered like eConnect™. The eCord is a standard feature on ACQUITY UPLC™ and some UHPLC Columns for use with Waters ACQUITY™ and Arc™ UHPLC Systems. The eConnect Column tag is an optional feature for HPLC columns and are designed to be used with the Waters Alliance™ iS HPLC System.

Although the eCord Chip and eConnect Electronic Column Tag designs differ, they are designed to be robust and durable. When in operation, the eCord enabled column is connected via a cable tether to the system, allowing the electronic tag to be physically located outside the column manager. This placement reduces exposure to variable operating temperatures and solvents. In contrast, the eConnect Tag is designed to be within the column manager, it is exposed to the full range of operating temperatures and may be exposed to solvents. Extensive testing has been conducted to ensure the performance of the eConnect Tag in this environment. The durability of the eConnect Tag was confirmed with various solvents and solutions with pH between 2 and 11. In the case of significant solvent exposure, communication between the tag and system can be temporarily interrupted. Allowing the eConnect tag to dry following solvent exposure, evaporation should be suitable in most cases, will allow the tag performance to recover.

In addition, eConnect Tags also offer a convenient, easy-to-use connection to the Waters website. Using an NFC enabled internet connected device to scan the eConnect Tag will provide the user with a web address which will take them directly to the Waters product page for that consumable. This page provides access to product information, like Care and Use or Certificate of Analysis, as well as reorder and related product support.

b. Use of Narrow-bore (<3.0 mm i.d.) Columns

This section describes how to minimize extra column effects and provides guidelines on maximizing the performance of a narrow-bore column in an HPLC system. A 3.0 mm i.d. narrow-bore column usually requires no system modifications. A 2.1 mm i.d. column, however, requires modifications to the HPLC system in order to eliminate excessive system bandspreading volume. Without proper system modifications, excessive system bandspreading volume causes peak broadening and has a large impact on peak width as peak volume decreases.

c. Impact of Bandspreading Volume on 2.1 mm i.d. Column Performance

Note: Flow splitters after the column will introduce additional bandspreading.

System optimization, especially in a system that contains a flow splitter, can have dramatic effects on sensitivity and resolution. Optimization includes using correct-depth ferrules and minimizing tubing diameter and lengths. An example is given in Figure 10 where system optimization resulted in a doubling of sensitivity and resolution of the metabolite in an LC/MS/MS system.

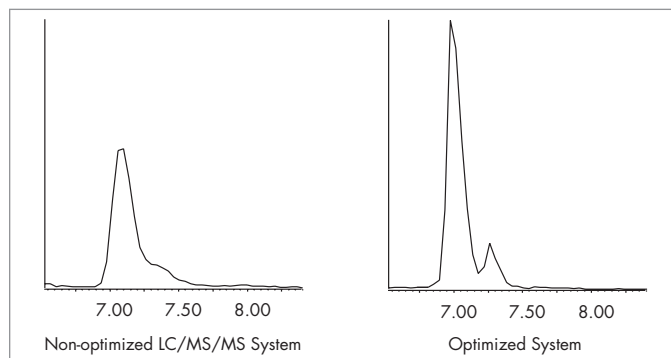


Figure 10: Non-Optimized vs. Optimized LC/MS/MS System.



d. Non-Optimized vs. Optimized LC/MS/MS System: System Modification Recommendations

1. Use a microbore detector flow cell with ≤ 2.1 mm i.d. columns.

Note: Detector sensitivity is reduced with the shorter flow cell path length in order to achieve lower bandspreading volume.

2. Minimize injector sample loop volume.
3. Use 0.009 inch (0.25 mm) tubing between pump and injector.
4. Use 0.009 inch (0.25 mm) tubing for rest of connections in standard systems and 0.005 inch (0.12 mm) tubing for narrowbore (2.1 mm i.d.) systems.
5. Use perfect (pre-cut) connections (with a variable depth inlet if using columns from different suppliers).
6. Detector time constants should be shortened to less than 0.2 seconds.

Waters™

Waters, ACQUITY, ACQUITY UPLC, eCord, eConnect, Empower, Alliance, Arc, SunFire, Oasis, and Sep-Pak are trademarks of Waters Corporation. Acrodisc is a trademark of Pall Corporation. PEEK is a registered trademark of Agilent Technologies. SLIPFREE is a trademark of Thermo Fisher Scientific Inc. All other trademarks are the property of their respective owners.

Waters Corporation
34 Maple Street
Milford, MA 01757 U.S.A.
T: 1 508 478 2000
F: 1 508 872 1990
[waters.com](https://www.waters.com)